INSTRUCTION MANUAL

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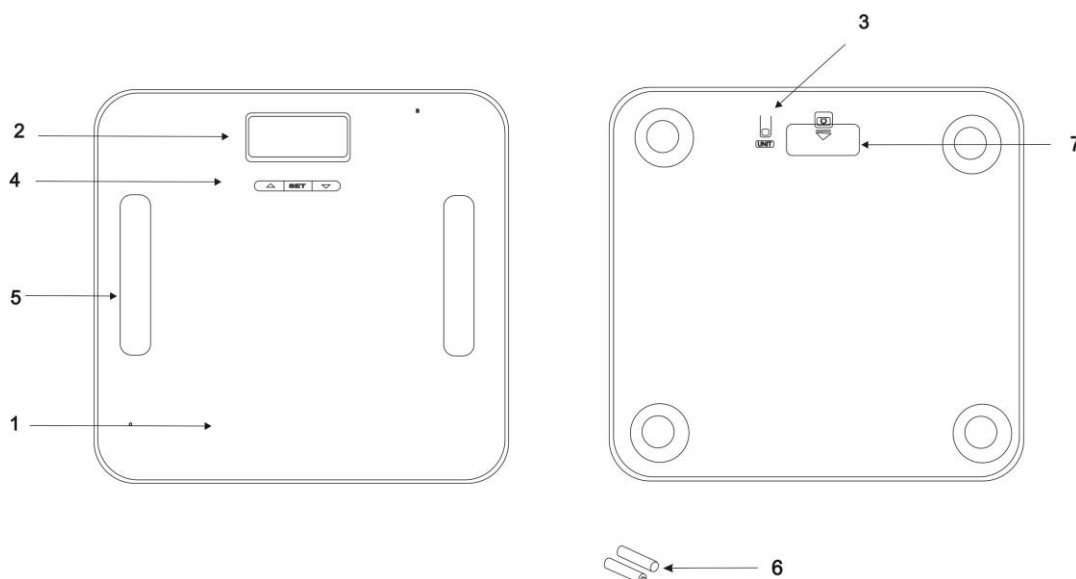
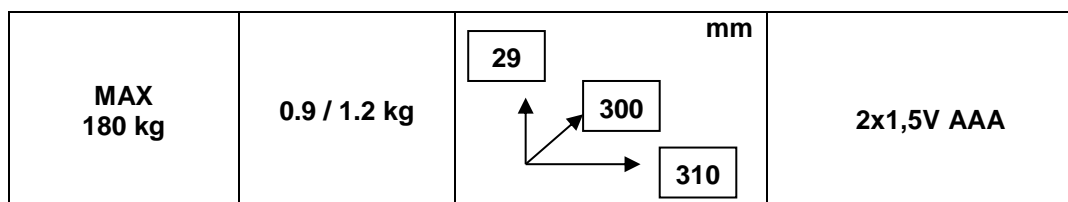
Indicators	Meaning
8888	Selection of the user and personal parameters
0000	Analysis of readings after weighing.
----	- Improper positioning of the scale. Place the scale on a flat, firm surface. - You are standing on the scale with your footwear on. Take off your footwear and step on the scale again.
LO	Indication of the low battery charge. The batteries are discharged, please, replace them with new ones.
ErrH	Indication of the high % of body fat. Body fat percentage is too high. Recommendations: you should keep to a low-calorie diet, do more exercise and pay more attention to your health.
Errl	Indication of the low % of body fat. Body fat percentage is too low. Recommendations: you should keep to a more nourishing diet and pay more attention to your health.
EEEE	Indication of overload. The weight exceeds maximum permissible value for measurement. Please, step off the platform to avoid damage to the scale.

CLEANING AND MAINTENANCE

- Wipe the scale with a soft cloth using a detergent and wipe dry.
- Do not use organic solvents, aggressive chemicals or abrasives.

STORAGE

- Observe all requirements from the CLEANING AND MAINTENANCE Section.
- Make sure there are no objects on the scale during storage.
- Store the scale in a cool dry place.

- The symbol on the unit, packing materials and/or documentations means used electrical and electronic units and battery's should not be tossed in the garbage with ordinary household garbage. These units should be passed to special receiving point.
- For additional information about actual system of the garbage collection address to the local authority.
- Valid utilization will help to save valuable resources and avoid negative work on the public health and environment which happens with incorrect using garbage.

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